

ANA BIO ISP

MULTICLEAN SOLUTION

Introduction

ANA Bio ISP Multiclean solution must be used on board of the class MIURA automatic analyzers for reaction cuvette washing.

Reagent storage & stability

The kit should be stored at the temperature +2-25°C. The expiry date of each reagent is indicated on its relative label. When stored at +2-25°C and protected from light. The solutions are stable until the expiry date stated on the label.

Preparation of the solution

Open the vial very carefully, avoiding liquid leakages. The content of one single vial (100 ml) of the multiclean solution must be diluted in to 5L of bi-distilled water in the tank for cleaner solution.

It is necessary that the solution, at the moment of the use, is at ambient temperature.

Operational procedure

See the instrument user manual.

Necessary material not furnished with the kit

- Automatic analyzer and/or for clinical chemistry class MIURA
- User manual of the relative instrument.

Warnings:

To obtain correct and reproducible results on the instrument, it is necessary to observe the following recommendations:








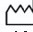





- Do not use solutions after their expiry date
- Avoid accurately their contamination.
- Avoid direct exposure to sunlight.

Precautions


To avoid personal and environmental contamination it is necessary to observe the following safety directives:

- Use protection glasses and disposable gloves during the manipulation
- Do not pipette reagent using your mouth
- Do not smoke, eat, drink or use cosmetics.
- The waste produced by the use of the product must be considered as dangerous waste and the dismantling of the same should be done according the local laws

Symbols

 IVD	In Vitro Diagnostics		Caution
 LOT	Batch No.		Product Expiry Date
 CONT	Content		Manufactured By
	Read Instructions		Date of Manufacture
	Storage Temperature		Keep Dry
 REF	Catalogue No.		Fragile
			Keep away from sun light



 KEE DIAGNOSTICS PVT LTD
(Formerly known as Kee GAD Biogen Pvt. Ltd.) .
CIN: U24231DL2004PTC128343 .
A subsidiary of KEE PHARMA LTD
A-8 . 3RD FLOOR. NARAINA INDUSTRIAL
AREA, PHASE 2 . NEW DELHI 110028.
W. www.keediagnosics.in.